

TECHNICAL DATA SHEET

CyFlow™ CD38 PE-DyLight™ 594 Anti-Hu; Clone HIT2

REF BN321669



Sysmex Partec GmbH

Arndtstraße 11 a-b

02826 Görlitz

Tel +49 3581 8746 0

E-mail: info@sysmex-partec.com

www.sysmex-partec.com

For Research Use Only.

Not for use in diagnostic or therapeutic procedures.

Specifications

Antigen	CD38
Alternative Names	T10
Clone	HIT2
Clonality	monoclonal
Format	PE-DyLight™ 594
Host / Isotype	Mouse / IgG1
Species Reactivity	Human
Negative Species Reactivity	—
Quantity	100 tests
Immunogen	Human thymocytes in foetus

The Safety Data Sheet for this product is available at www.sysmex-partec.com/services



Specificity

The mouse monoclonal antibody HIT2 recognizes CD38 antigen, a 45 kDa type II transmembrane glycoprotein strongly expressed mainly on plasma cells and activated T and B lymphocytes; it is an antigenic marker of lymphoid cells.

HLDA III—WS Code T 155

Application

The reagent is designed for Flow Cytometry analysis of human blood cells. Recommended usage is 4- μ l reagent-/ 100- μ l of whole blood or 10^6 cells in a suspension. The content of a vial (0.4 ml) is sufficient for 100 tests.

Other usages may be determined from the scientific literature.

Storage Buffer

The reagent is provided in stabilizing phosphate buffered saline (PBS) solution, pH \approx 7.4, containing 0.09% (w/v) sodium azide.

Storage and Stability

Storage	Avoid prolonged exposure to light. Store in the dark at 2-8°C. Do not freeze.
Stability	Do not use after expiration date stamped on vial label.

Background Information

CD38 (NAD⁺ glycohydrolase) is a type II transmembrane glycoprotein able to induce activation, proliferation and differentiation of mature lymphocytes and mediate apoptosis of myeloid and lymphoid progenitor cells. Another role of CD38 is provided by enzymatic activity of its extracellular part. CD38 acts as NAD⁺ glycohydrolase converting NAD⁺ into ADP-ribose, as ADP-ribosyl cyclase producing cADPR and as cADPR hydrolase, thus affecting levels of calcium-mobilizing metabolites. ADPR produced by CD38 serves as an important second messenger of neutrophil and dendritic cell migration.

Warnings

Non-Hazardous Statement: This is not considered hazardous by the criteria in 29 CFR 1910.1200 or the General Classification guideline for preparations of the EU.




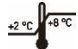





Safety Data Sheet Statement: Important information regarding the safe handling, transport, and disposal of this product is contained in the Safety Data Sheet.

The Safety Data Sheet for this product is available at www.sysmex-partec.com/services

References

- McMichael AJ, Beverley PCL, Cobbold S, et al (Eds): Leucocyte Typing III, White Cell Differentiation Antigens. Oxford University Press, Oxford. 1987; 1-1050. < NLM ID: 8913266 >
- Cakir-Kiefer C, Muller-Steffner H, Oppenheimer N, Schuber F: Kinetic competence of the cADP-ribose-CD38 complex as an intermediate in the CD38/NAD⁺ glycohydrolase-catalysed reactions: implication for CD38 signalling. *Biochem J.* 2001 Sep 1; 358(2):399-406. < PMID: 11513738 >
- Lund FE, Muller-Steffner H, Romero-Ramirez H, Moreno-García ME, Partida-Sánchez S, Makris M, Oppenheimer NJ, Santos-Argumedo L, Schuber F: CD38 induces apoptosis of a murine pro-B leukemic cell line by a tyrosine kinase-dependent but ADP-ribosyl cyclase- and NAD glycohydrolase-independent mechanism. *Int Immunol.* 2006 Jul; 18(7):1029-42. < PMID: 16720618 >
- Kolar GR, Mehta D, Pelayo R, Capra JD: A novel human B cell subpopulation representing the initial germinal center population to express AID. *Blood.* 2007 Mar 15; 109(6):2545-52. < PMID: 17132718 >
- Partida-Sanchez S, Gasser A, Fliegert R, Siebrands CC, Dammermann W, Shi G, Mousseau BJ, Sumoza-Toledo A, Bhagat H, Walseth TF, Guse AH, Lund FE: Chemotaxis of mouse bone marrow neutrophils and dendritic cells is controlled by adp-ribose, the major product generated by the CD38 enzymereaction. *J Immunol.* 2007 Dec 1; 179(11):7827-39. < PMID: 18025229 >
- Rozková D, Novotná L, Pytlík R, Hochová I, Kozák T, Bartůnková J, Spísek R: Toll-like receptors on B-CLL cells: expression and functional consequences of their stimulation. *Int J Cancer.* 2010 Mar 1; 126(5):1132-43. < PMID: 19685493 >

Symbols

	Reagent catalogue number		Contains sufficient for <n> tests
	For Research Use Only. Not for use in diagnostic and therapeutic procedures.		Temperature limitation
	Batch code		Keep away from sunlight
	Manufacturer		Consult accompanying documents
	Use by		

This product is provided under an intellectual property license from Life Technologies Corporation. The transfer of this product is conditioned on the buyer using the purchased product solely in research conducted by the buyer, excluding contract research or any fee for service research, and the buyer must not sell or otherwise transfer this product or its components for (a) diagnostic, therapeutic or prophylactic purposes; (b) testing, analysis or screening services, or information in return for compensation on a per-test basis; (c) manufacturing or quality assurance or quality control, or (d) resale, whether or not resold for use in research. For information on purchasing a license to this product for purposes other than as described above, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.